




AWC®duo – Assembly Instructions


	Ref. № 200.57.05	Ref. № 200.57.06
Endoscope diameter	9.8–11.5 mm	11.5–13.5 mm
Insertion length	103–110 cm	160–170 cm

Assembly
video

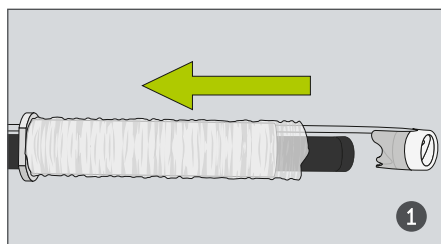


Slide the endoscope through the small tube supplied for the endoscope sleeve.


Important: moisten the endoscope tip and the flexible cap with water.

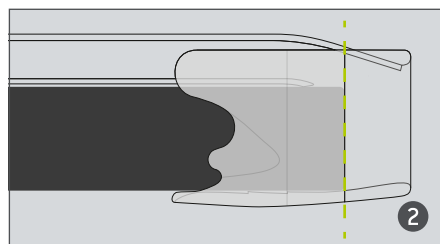


Slide the flexible cap on the endoscope tip using the strap until the endoscope tip reaches the transparent part of the cap.

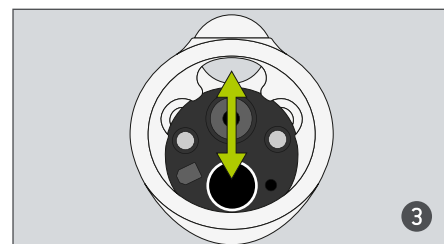


Position the AWC®duo working channel suitably for the application. For maximum distance, place the working channels on opposite sides.

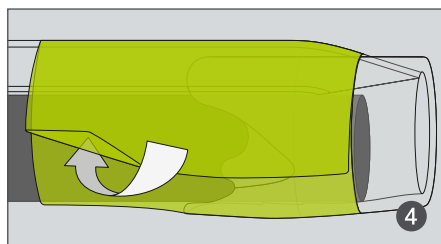


Fasten a strip of adhesive tape on the translucent part of the cap and the endoscope shaft. Fold one edge for easy dismantling.

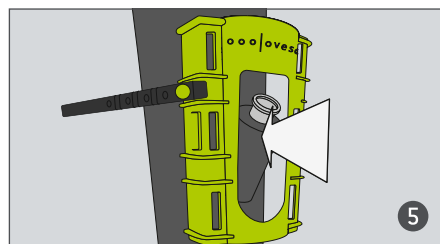


Remove the inlet valve from the endoscope working channel. Push the mounting clasp onto the endoscope.

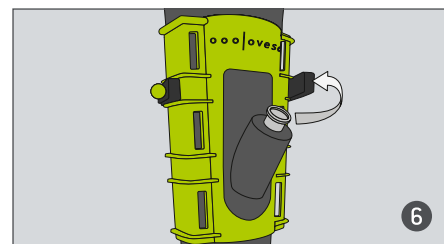


Fix the mounting clasp with the black fastening strap.

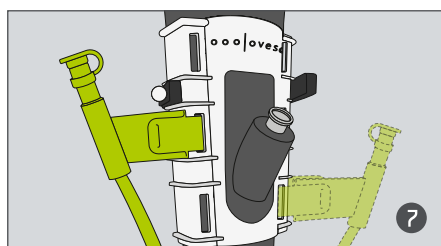


Fasten the connector in a suitable position, either on the left or right side, depending on the user's preference. Close the inlet valve of the additional working channel.

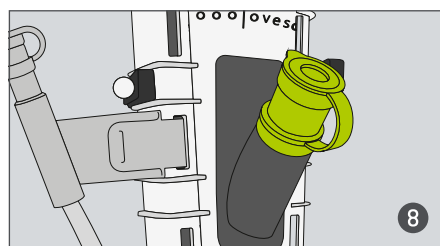


Re-mount the closed inlet valve in the working channel entrance of the endoscope.

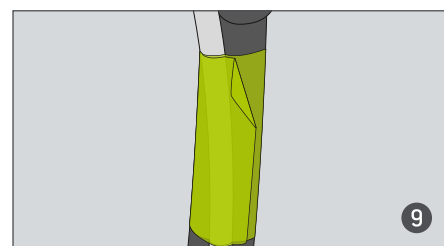


Attach the AWC®duo working channel below the endoscope handle to the endoscope shaft using one strip of adhesive tape.

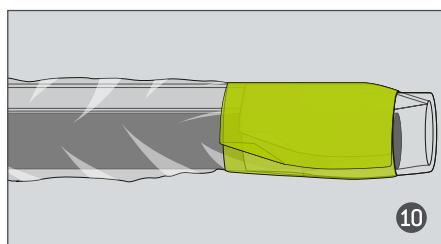

Pull the sleeve end to the flexible cap and position it at the strap. Place adhesive tape over the existing tape to secure the sleeve to the flexible cap.

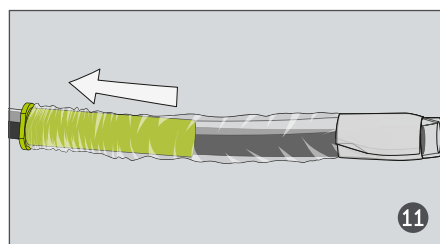

Pull the endoscope sleeve back to the proximal end of the endoscope using the small tube. Completely unroll the sleeve from the small tube.

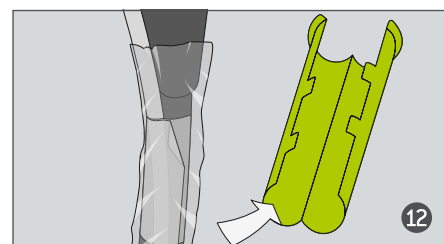

Remove the small tube. Smooth out the sleeve towards the proximal end to remove excess material at the distal end.