


ovesco training campus
TÜBINGEN
CARY
PARIS
CHENNAI


OTSC®neo System – Assembly instructions

Endoscope diameter			
OTSC®neo mini	OTSC®neo 11	OTSC®neo 12	OTSC®neo 14
8.5 – 10 mm	8.5 – 11 mm	10.5 – 12 mm	11.5 – 14 mm


Assembly video

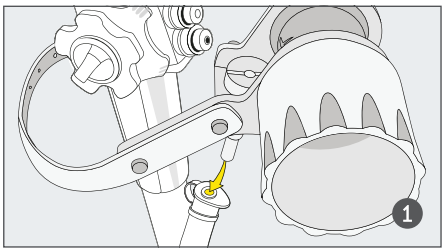


Mount hand wheel on the working channel of the endoscope.

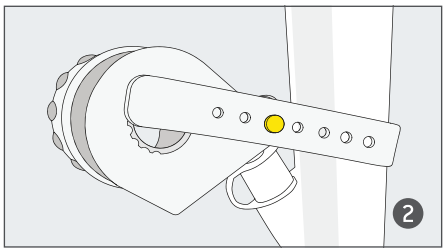


Secure the fastening strap of the hand wheel in the corresponding hole.

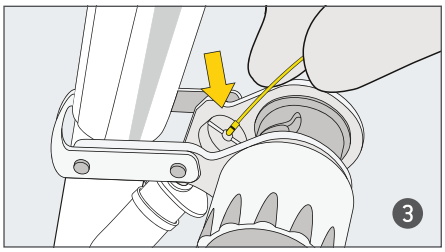


Insert thread retriever into the working channel.

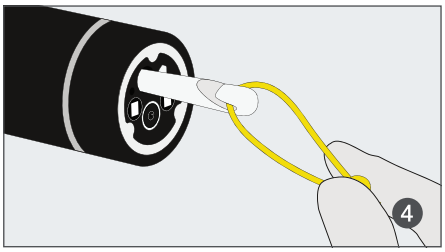


Grasp the end of the thread with the thread retriever. Pull the thread through the working channel.

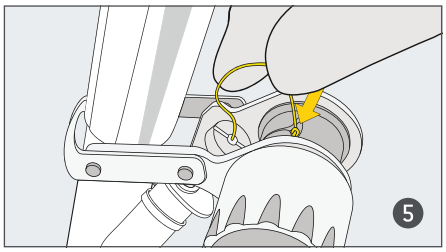


Clamp the thread in the hand wheel's corresponding gap.

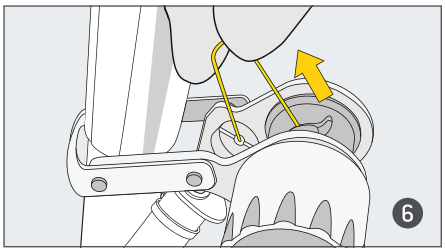


Hold and pull the thread.

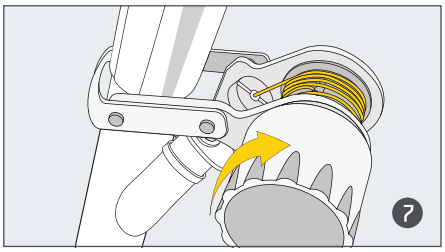


Wind the thread with the hand wheel.

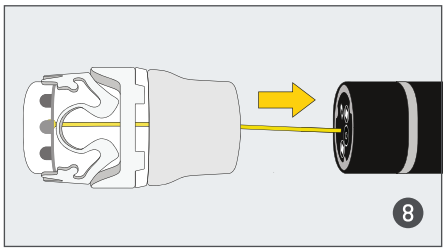


Place cap onto the endoscope.

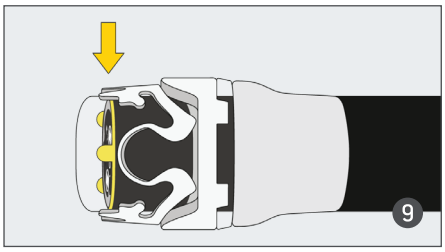


Push endoscope forward until it abuts against the stoppers in the cap.
Aids: fluids, e.g. gel.

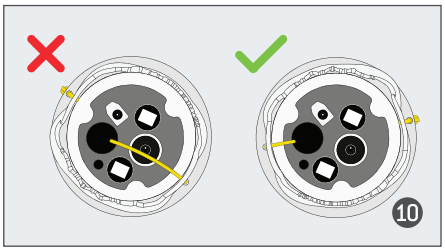

Mind position of the thread while attaching the cap onto the endoscope. The thread must not interfere with the field of vision.

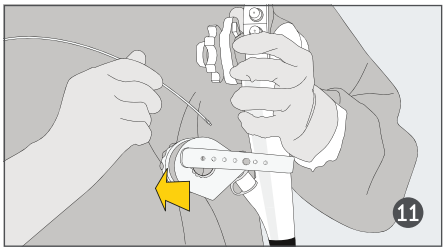

If another instrument is inserted into the working channel, the thread needs to be mildly tensioned.

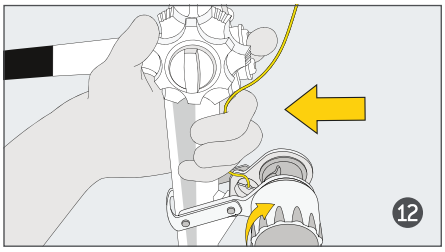

Before clipping, the instrument needs to be secured. To release the clip, the hand wheel needs to be twisted quickly and completely.

OVE_OTSCneo-System_infosheet-assembly_eng_Rev02

Ovesco Endoscopy AG | Friedrich-Miescher-Str. 9 | 72076 Tuebingen | Germany | +49 (0)7071 96528-160 | service@ovesco.com | www.ovesco.com